



activ5™


Genuine Accessories

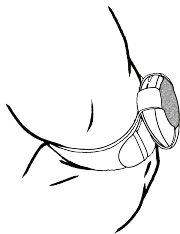

USER MANUAL

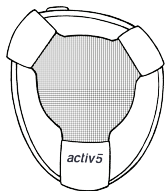


Please use ActivBody products and accessories solely for the designated purposes and methods of use, after consulting with your responsible physician, therapist or trainer and following the directions on the associated manuals and software applications.

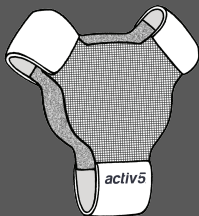


Accessory Attachment


The Accessory Attachment rigs each side of the Activ5 device with a contact surface disk (either made of Hooks or Loops) with which other compatible Loops or Hooks surfaced accessories or objects can be attached.

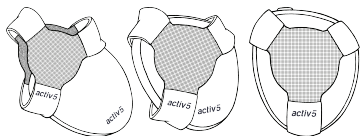

- For proper fit, slip the Activ5 device into the Accessory Attachment so that the Activ5 logos on both items match their positions.

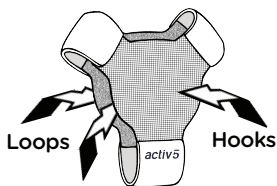


- If the Accessory Attachment's logo is visible (facing out), there will be the Hooks disk next to it and a Loops disk on the opposite side of the Activ5. This configuration is useful for attaching the Activ5 to the Wrap or other Loops surfaced accessories or objects.

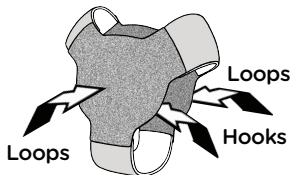

- Flipping the Accessory Attachment inside out provides comfort and more combinations for attaching the Activ5 with Hooks surfaced accessories or objects.



Do not stretch the Accessory Attachment more than just enough to sneak in and adjust the Activ5 Device!

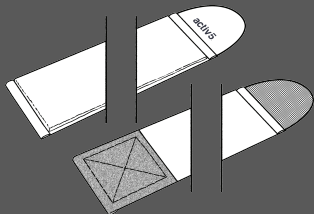

Wrap

The Wrap enables attaching the Activ5 device to different body parts - like palms, wrists, elbows ankles, feet, etc.

